


Certificate of Analysis
For Laboratory Use Only

PRODUCT NAME	Ursolic acid
PART NUMBER	00021600
STANDARD TYPE	Primary (P)
LOT NUMBER	00021600-022
REPORT NUMBER	CDXA-RSS-3985-03
SAMPLE NUMBER	CDXA-10-2823
DATE OF SAMPLE	05/05/2010
DATE OF RE-EVALUATION	06/09/2014 (1 st); 08/15/2018 (2 nd); 07/07/2022 (3 rd)
DATE OF REPORT	08/22/2022

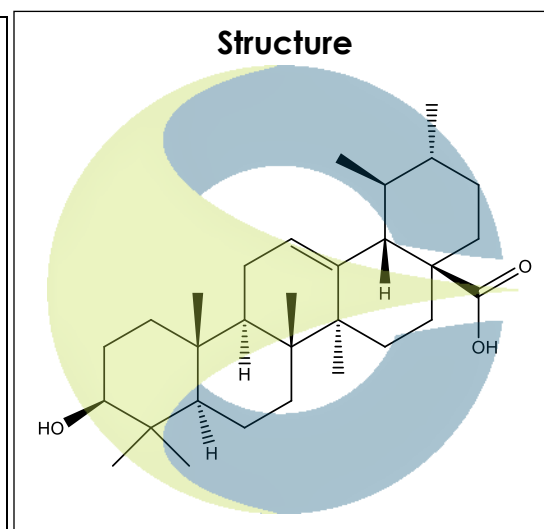

CHEMICAL NAMES	(3β)-3-Hydroxyurs-12-en-28-oic acid; Bungeolic acid; Formosolic acid; Malol; Prunol; Urson
CHEMICAL FORMULA	C ₃₀ H ₄₈ O ₃
MOLECULAR WEIGHT (MW)	456.70
CHEMICAL FAMILY	Triterpenes
CAS NUMBER	[77-52-1]
EC#(EINECS)	201-034-0
RTECS	YU9520000

ANALYTICAL RESULTS

TEST	METHOD	SPECIFICATION	RESULT
Adjusted Purity	NA	≥ 85%	95.1%
HPLC	0.700.10.2.METH22	≥ 90%	96.7%
NMR	0.700.12.4	Conforms to structure	Conforms to structure
Mass Spectrum	Direct Infusion, ESI (-)	Conforms	Conforms
Residual Solvent*	99.1-CD-3.0-000115	NA	EtOH -1.2%
Water	0.700.12.37	NA	0.5%
Appearance	NA	NA	White powder

*Testing performed at an ISO 17025 accredited subcontracted laboratory.

ADJUSTED PURITY: 95.1% IS BASED ON (100 - 1.2 SOLVENTS - 0.5 WATER) X 96.7% HPLC **STORAGE CONDITIONS**

STORAGE Room Temperature in a dry place.

EXPIRATION DATE 07/2027 under the above conditions.

This document is the property of ChromaDex, Inc. and/or its relevant affiliates and contains confidential and proprietary material for the sole use of the intended recipient(s). Any review, use, distribution or disclosure by or to others is strictly prohibited.

10900 Wilshire Blvd., Suite 600, Los Angeles, CA 90024 | T: +1 949-419-0288 | F: +1 949-419-0294 | www.chromadex.com
 Laboratory address: 1751 S. Fordham St., Suite 350, Longmont, CO 80503

India Contact: Life Technologies (India) Pvt. Ltd.
 306, Aggarwal City Mall, Opposite M2K Pitampura, Delhi - 110034 (INDIA).
 Ph: +91-11-42208000, 42208111, 42208222, Mobile: +91-9810521400, Fax: +91-11-42208444
 Email: customerservice@lifetechindia.com Website: www.lifetechindia.com


Certificate of Analysis

For Laboratory Use Only

ANALYTICAL CONDITIONS

INSTRUMENT	AGILENT 1260 HPLC UV-VIS (DAD) CONNECTED TO DIONEX CORONA VEO CHARGED AEROSOL DETECTOR (CAD); AGILENT 1100 ION TRAP
COLUMN	Phenomenex Luna C18(2) 150 x 2.00 mm, 5 µm particle size
MOBILE PHASE	A - 0.1% Trifluoroacetic acid in Ultrapure water, B - Acetonitrile; Isocratic 5% B for 5 minutes, then increasing to 95% B over 20 minutes. Hold at 95% B for 15 minutes, then decreasing to 5% B over 1 minute.
COLUMN TEMP.	60 °C
FLOW RATE	0.4 mL/minute
INJECTION VOL.	2.0 µL
INJECTION CONC.	1.0 mg/mL in Methanol
DETECTION	CAD

HPLC CHROMATOGRAM OF URSOLIC ACID (CDXA-10-2823)

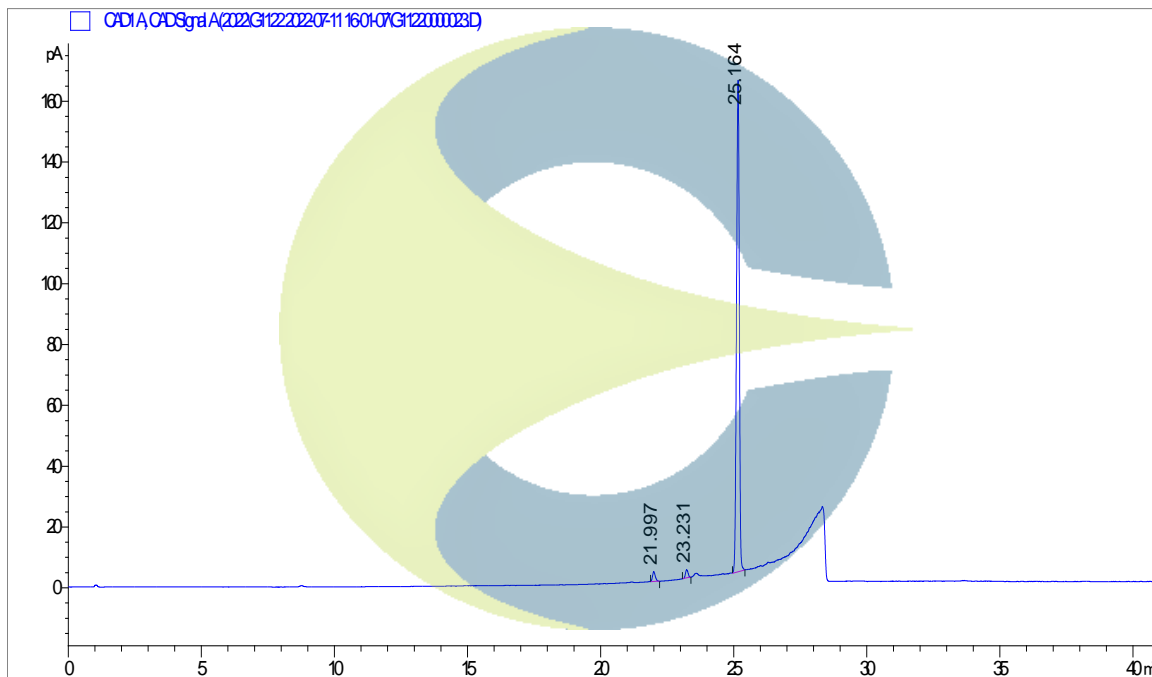


This document is the property of ChromaDex, Inc. and/or its relevant affiliates and contains confidential and proprietary material for the sole use of the intended recipient(s). Any review, use, distribution or disclosure by or to others is strictly prohibited.


10900 Wilshire Blvd., Suite 600, Los Angeles, CA 90024 | T: +1 949-419-0288 | F: +1 949-419-0294 | www.chromadex.com
 Laboratory address: 1751 S. Fordham St., Suite 350, Longmont, CO 80503

India Contact: Life Technologies (India) Pvt. Ltd.
 306, Aggarwal City Mall, Opposite M2K Pitampura, Delhi - 110034 (INDIA).
 Ph: +91-11-42208000, 42208111, 42208222, Mobile: +91-9810521400, Fax: +91-11-42208444
 Email: customerservice@lifetechindia.com Website: www.lifetechindia.com

HPLC CHROMATOGRAM OF DILUENT BLANK

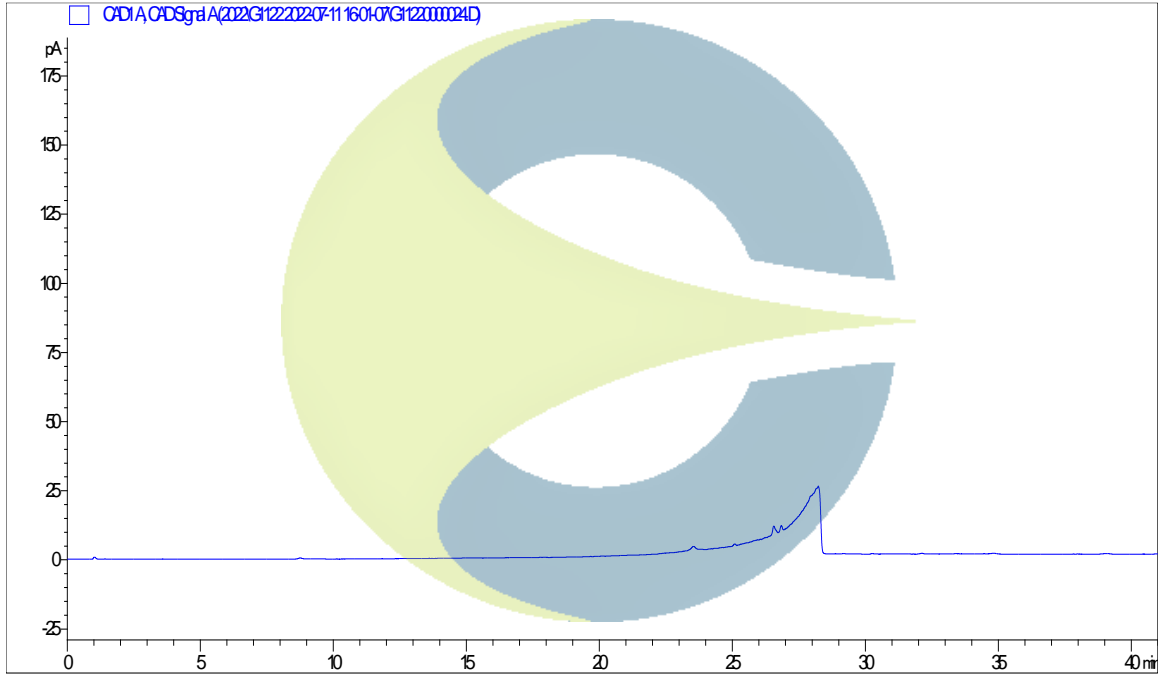

MASS SPECTRUM OF URSOLIC ACID (CDXA-10-2823)

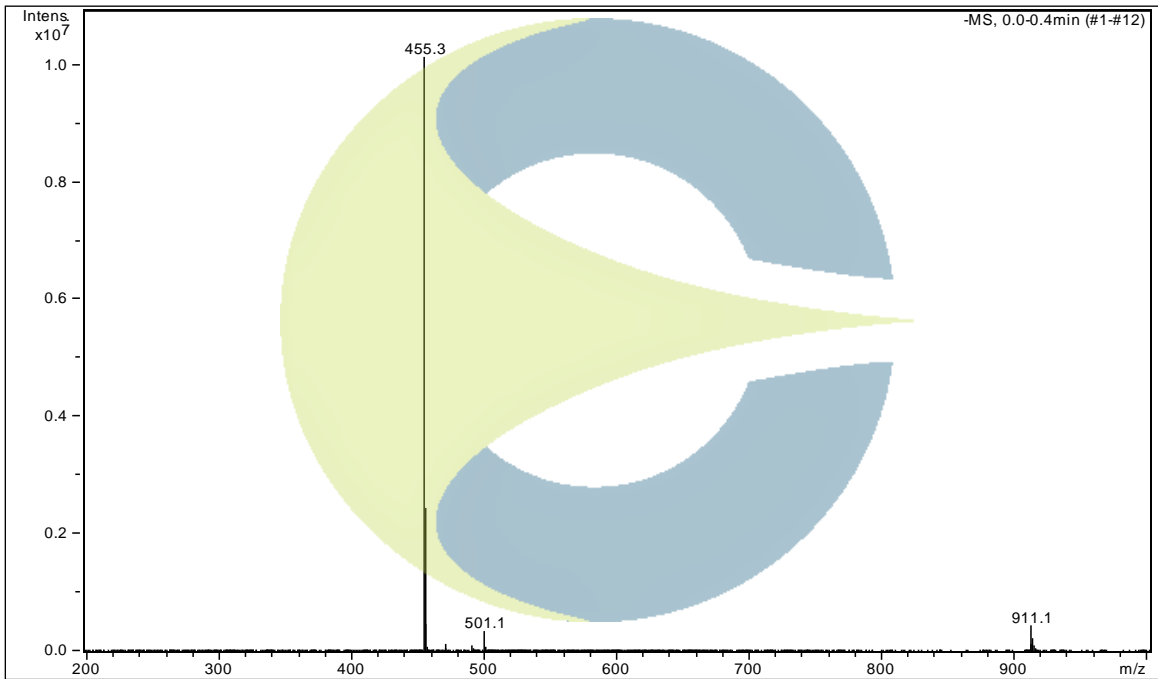


This document is the property of ChromaDex, Inc. and/or its relevant affiliates and contains confidential and proprietary material for the sole use of the intended recipient(s). Any review, use, distribution or disclosure by or to others is strictly prohibited.


10900 Wilshire Blvd., Suite 600, Los Angeles, CA 90024 | T: +1 949-419-0288 | F: +1 949-419-0294 | www.chromadex.com
Laboratory address: 1751 S. Fordham St., Suite 350, Longmont, CO 80503

India Contact: Life Technologies (India) Pvt. Ltd.
306, Aggarwal City Mall, Opposite M2K Pitampura, Delhi - 110034 (INDIA).
Ph: +91-11-42208000, 42208111, 42208222, Mobile: +91-9810521400, Fax: +91-11-42208444
Email: customerservice@lifetechindia.com Website: www.lifetechindia.com

¹H-NMR SPECTRUM OF URSOLIC ACID (CDXA-10-2823) IN PYRIDINE-*d*₅

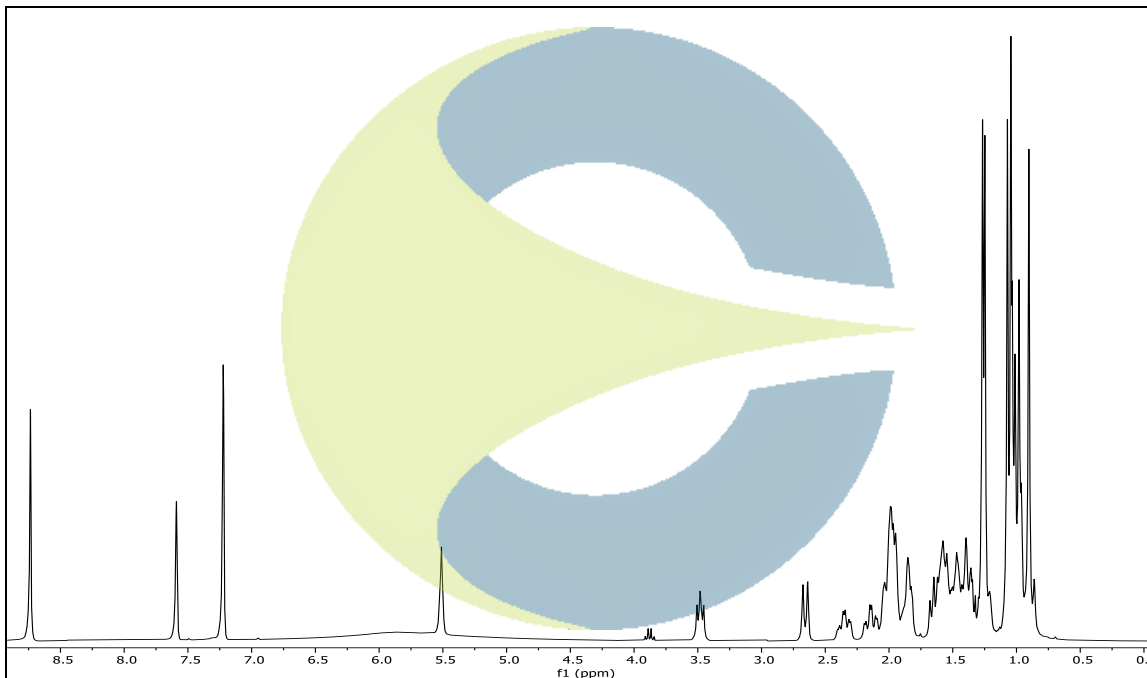


¹³C-NMR SPECTRUM OF URSOLIC ACID (CDXA-10-2823) IN PYRIDINE-*d*₅

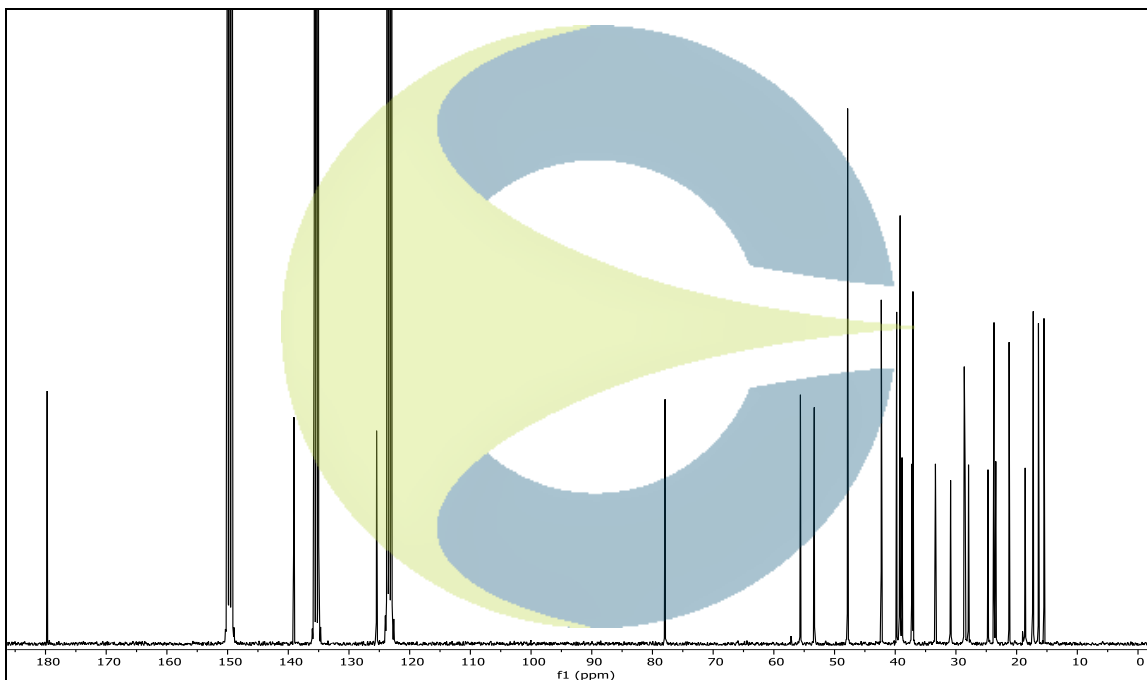

This document is the property of ChromaDex, Inc. and/or its relevant affiliates and contains confidential and proprietary material for the sole use of the intended recipient(s). Any review, use, distribution or disclosure by or to others is strictly prohibited.

10900 Wilshire Blvd., Suite 600, Los Angeles, CA 90024 | T: +1 949-419-0288 | F: +1 949-419-0294 | www.chromadex.com
Laboratory address: 1751 S. Fordham St., Suite 350, Longmont, CO 80503

India Contact: Life Technologies (India) Pvt. Ltd.
306, Aggarwal City Mall, Opposite M2K Pitampura, Delhi - 110034 (INDIA).
Ph: +91-11-42208000, 42208111, 42208222, Mobile: +91-9810521400, Fax: +91-11-42208444
Email: customerservice@lifetechindia.com Website: www.lifetechindia.com

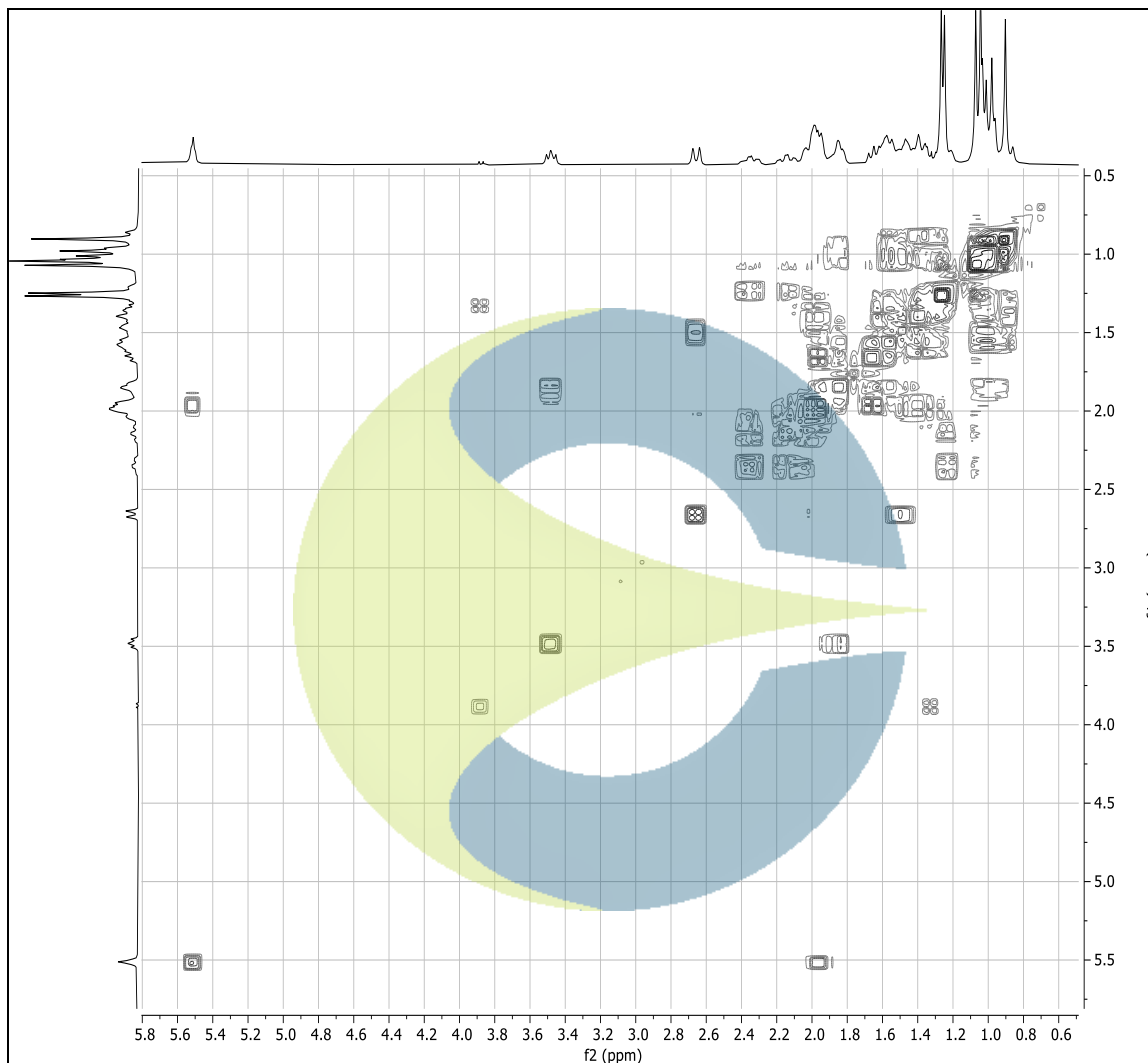COSY-NMR SPECTRUM OF URSOLIC ACID (CDXA-10-2823) IN PYRIDINE-*d*₅

REVISION HISTORY

<u>Revision History</u>	<u>Date of Revision</u>	<u>Document/Changes</u>
00	05/20/2010	New report
01	06/12/2014	Passed reevaluation; updated expiration date
02	09/22/2018	Passed second reevaluation. Updated analytical results, analytical conditions, chemical structure format, expiration date and chromatograms. Removed melting point.
03	08/22/2022	Passed third re-evaluation by HPLC, KF, NMR. Updated expiration date, analytical results, analytical conditions, and chromatograms.

This document is the property of ChromaDex, Inc. and/or its relevant affiliates and contains confidential and proprietary material for the sole use of the intended recipient(s). Any review, use, distribution or disclosure by or to others is strictly prohibited.

10900 Wilshire Blvd., Suite 600, Los Angeles, CA 90024 | T: +1 949-419-0288 | F: +1 949-419-0294 | www.chromadex.com
 Laboratory address: 1751 S. Fordham St., Suite 350, Longmont, CO 80503

India Contact: Life Technologies (India) Pvt. Ltd.
 306, Aggarwal City Mall, Opposite M2K Pitampura, Delhi - 110034 (INDIA).
 Ph: +91-11-42208000, 42208111, 42208222, Mobile: +91-9810521400, Fax: +91-11-42208444
 Email: customerservice@lifetechindia.com Website: www.lifetechindia.com