


Certificate of Analysis
For Laboratory Use Only

PRODUCT NAME	Cinnamic acid methylester
PART NUMBER	00003664
STANDARD TYPE	Analytical Standard (AS)
LOT NUMBER	00003664-13V
REPORT NUMBER	CDXA-RSS-4497-02
SAMPLE NUMBER	CDXA-11-2063
DATE OF SAMPLE	03/29/2011
DATE OF RE-EVALUATION	08/03/2015 (1 st); 06/09/2020 (2 nd)
DATE OF REPORT	07/06/2020

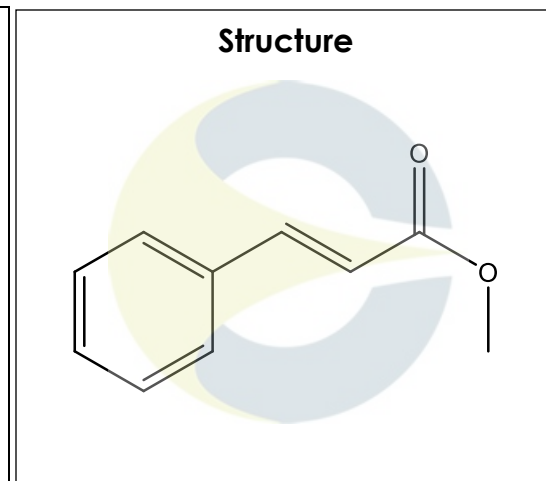

CHEMICAL NAMES	3-Phenyl-2-propenoic acid methyl ester; 3-Phenylacrylic acid methyl ester; Methyl 3-phenyl-2-propenoate; Methyl 3-phenylacrylate; Methyl cinnamate
CHEMICAL FORMULA	C ₁₀ H ₁₀ O ₂
MOLECULAR WEIGHT (MW)	162.19
CHEMICAL FAMILY	Phenolic acids
CAS NUMBER	[103-26-4]
EC#(EINECS)	203-093-8
RTECS	GE0190000

ANALYTICAL RESULTS

TEST	METHOD	SPECIFICATION	RESULT
GC	0.700.10.2.METH140	NA	99.9%
Mass Spectrum	0.700.10.2.METH140	Conforms	Conforms
Appearance	NA	NA	White powder


STORAGE CONDITIONS

STORAGE	Room Temperature in a dry place.
EXPIRATION DATE	06/2025 under the above conditions.

This document is the property of ChromaDex, Inc. and/or its relevant affiliates and contains confidential and proprietary material for the sole use of the intended recipient(s). Any review, use, distribution or disclosure by or to others is strictly prohibited.


ANALYTICAL CONDITIONS

INSTRUMENT	AGILENT 6890N WITH FLAME IONIZATION DETECTOR (FID) AND AGILENT 5973N MASS SELECTIVE DETECTOR (MSD)
COLUMN	Phenomenex Zebron ZB-5HT Inferno, 30 m x 0.32 mm x 0.25 µm; connected to an MXT Y-union, connected to a Restek Siltek deactivated guard column: 10 m x 0.25 mm to FID; and a 5 m x 100 µm transfer line to MSD.
CARRIER GAS	Helium
COLUMN TEMP.	Isothermic 40 °C for 5.0 minutes, then increasing at 10.0 °C/minute to 340 °C. Hold at 340 °C for 5.0 minutes.
FLOW RATE	2.0 mL/minute
INJECTOR TEMP.	300 °C
INJECTOR SPLIT RATIO	20:1
INJECTION VOL.	1.0 µL
INJECTION CONC.	1.5 mg/mL in Methanol
DETECTION	FID (350 °C) & MSD (Transfer line 300 °C, Source 250 °C, Quad 200 °C)

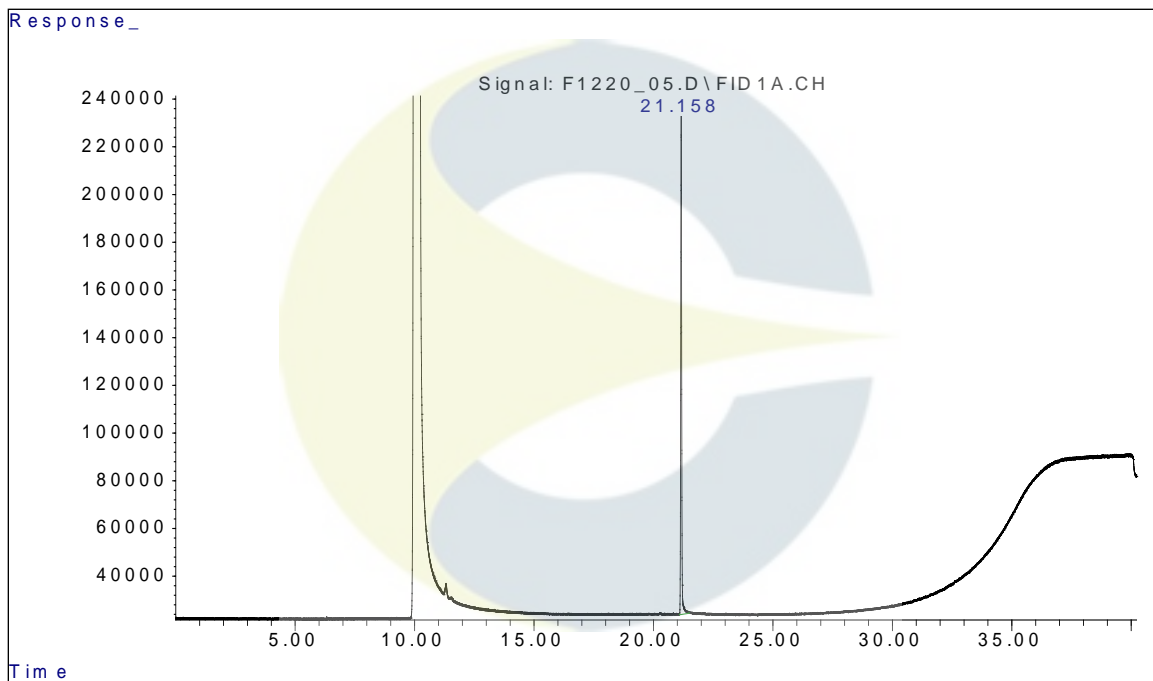
GC CHROMATOGRAM OF CINNAMIC ACID METHYLESTER (CDXA-11-2063)


This document is the property of ChromaDex, Inc. and/or its relevant affiliates and contains confidential and proprietary material for the sole use of the intended recipient(s). Any review, use, distribution or disclosure by or to others is strictly prohibited.

GC CHROMATOGRAM OF DILUENT BLANK

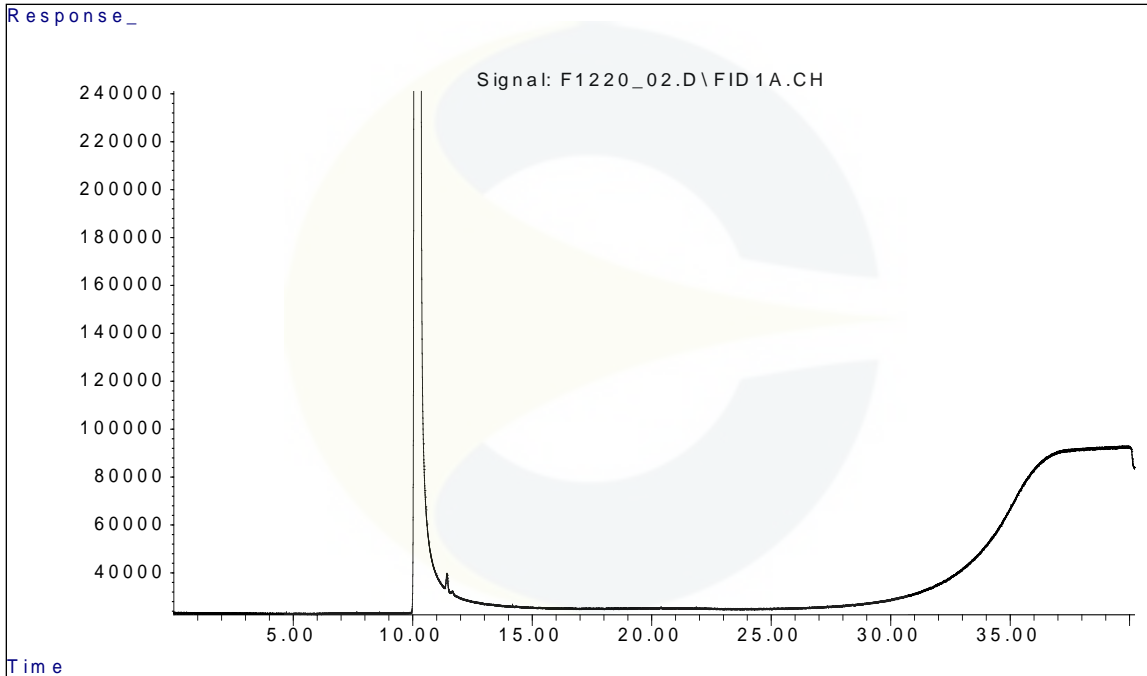

MASS SPECTRUM OF CINNAMIC ACID METHYLESTER (CDXA-11-2063)

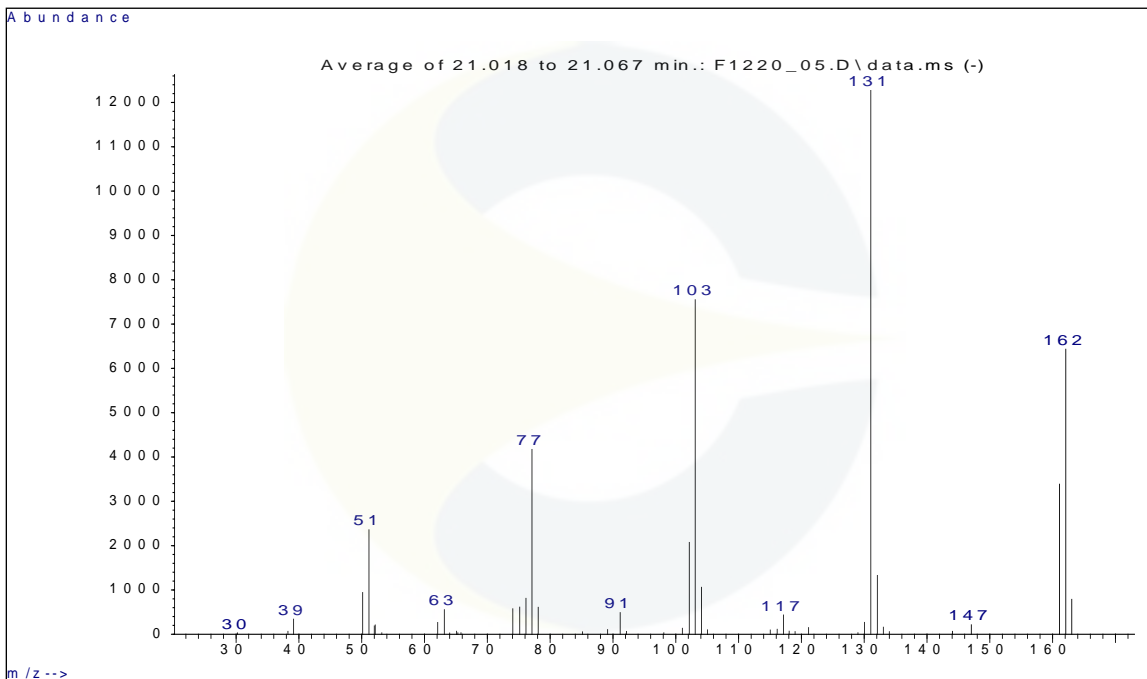

This document is the property of ChromaDex, Inc. and/or its relevant affiliates and contains confidential and proprietary material for the sole use of the intended recipient(s). Any review, use, distribution or disclosure by or to others is strictly prohibited.

REVISION HISTORY

<u>Revision History</u>	<u>Date of Revision</u>	<u>Document/Changes</u>
00	07/26/2011	New report
01	08/21/2015	Passed re-evaluation by GC/FID. Updated expiration date, standard type and melting point.
02	07/06/2020	Passed second re-evaluation by GC/FID MS. Updated expiration date, analytical method, conditions, chromatograms and results. Removed melting point.

This document is the property of ChromaDex, Inc. and/or its relevant affiliates and contains confidential and proprietary material for the sole use of the intended recipient(s). Any review, use, distribution or disclosure by or to others is strictly prohibited.