





SPECIMEN

For Laboratory Use Only

Certificate of Analysis

| PRODUCT | Morinda citrifolia (Noni) Fruit |
|----------------|------------------------------------|
| PART NUMBER | 00030846 |
| REFERENCE TYPE | Biomass Reference Material (BRM) * |
| LOT NUMBER | 00030846-421 |
| COMMON NAME | Noni |
| LATIN NAME | Morinda citrifolia L. [Rubiaceae] |
| PLANT PART | Fruit |
| SAMPLE NUMBER | CDXP-21-00003 |
| REPORT NUMBER | CDXA-BRM-198-00 |
| DATE OF SAMPLE | 01/05/2021 |
| DATE OF REPORT | 02/01/2021 |

^{*}Note - Biomass Reference Material (BRM) is not a voucher specimen

ANALYTICAL RESULTS

| TEST | METHOD | RESULT |
|--------------------|---|--|
| Appearance | Macroscopy | (1,2) Whole/Cut: fruit ovoid to subglobose,2.5-5 cm diameter(3) Milled: Brown powder |
| Anatomy/Morphology | Microscopy (400x - acidified chloral hydrate glycerol solution) | (1) acicular calcium oxalate crystals |
| HP-TLC | See HP-TLC conditions | Conforms |

STORAGE CONDITIONS

STORAGE 20-30 °C; Dry storage area; Insect free; Volatile free

EXPIRATION DATE 01/2026 under the above conditions

This document is the property of ChromaDex, Inc. and/or its relevant affiliates and contains confidential and proprietary material for the sole use of the intended recipient(s). Any review, use, distribution, or disclosure by or to others is strictly prohibited.







SPECIMEN

For Laboratory Use Only

Certificate of Analysis

MACROSCOPY

(1)

1)






(2)

MICROSCOPY

(3)



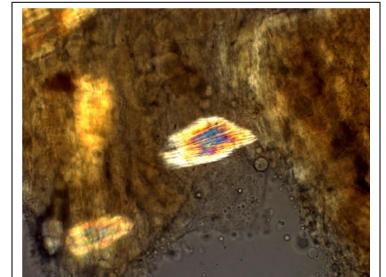

(1)

This document is the property of ChromaDex, Inc. and/or its relevant affiliates and contains confidential and proprietary material for the sole use of the intended recipient(s). Any review, use, distribution, or disclosure by or to others is strictly prohibited.







SPECIMEN

Certificate of Analysis

For Laboratory Use Only

HP-TLC CONDITIONS

STATIONARY PHASE Silica gel 60, HPTLC plates

SAMPLE PREPARATION 0.3 g + 3 mL Methanol, sonicate/heat at 50 °C for 30 min.

Chloroform/ Methanol/ Formic acid/ Water **MOBILE PHASE**


[6/2/1.4/0.1]

DETECTION (1) 10% Sulfuric, 100°C, 2min, Visible light


(2) 10% Sulfuric, 100°C, 2min, 366nm

HP-TLC PLATES



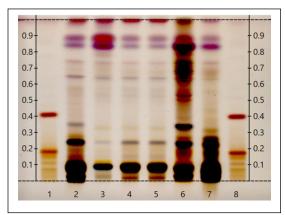

HP-TLC (2)

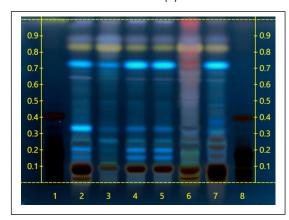

ID

HP-TLC LANE APPLICATIONS

| Lane | ID | Lane | ID |
|------|-----------------------------------|------|-------------------------------------|
| 1 | Epigallocatechin (3 µL) | 5 | Morinda citrifolia fruit BRM (3 µL) |
| 2 | Morinda citrifolia (fruit) (3 μL) | 6 | Morinda citrifolia (leaf) (3 μL) |
| 3 | Morinda citrifolia (fruit) (3 μL) | 7 | Morinda citrifolia (pulp) (3 μL) |
| 4 | Morinda citrifolia BRM (3 μL) | 8 | Epigallocatechin (3 µL) |

REVISION HISTORY

| <u>Revision History</u> | <u>Date of Revision</u> | <u>Document/Changes</u> |
|-------------------------|-------------------------|-------------------------|
| | | |

00 02/01/2021 New report

This document is the property of ChromaDex, Inc. and/or its relevant affiliates and contains confidential and proprietary material for the sole use of the intended recipient(s). Any review, use, distribution, or disclosure by or to others is strictly prohibited.